



Metro Center for High Technology Bldg. 2727 Second Ave. Suite 4113

Detroit, MI 48201

Phone: (313) 961-1606; Fax: (313)963-7130

Email: info@DetroitRandD.com
Web: www.DetroitRandD.com


## Prostaglandin E2 (PGE2) ELISA Kit with Immunoaffinity Columns

Cat # PGE1-IA: ELISA kit for measuring 8-isoprostane in biological and environmental samples: \$497

This is the same ELISA kit as our Cat# PGE1 kit only we have added 10 immunoaffinity columns for one-step purification of PGE2. Each immunoaffinity column can bind up to 1 microgram of PGE2. No extraction with organic solvents is necessary.

Buy in Quantity and Save! 2-4 kits \$472 each; 5-9 kits \$447 each; 10 or more kits \$422 each

This competitive ELISA kit is for determination of  $PGE_2$  (Prostaglandin  $E_2$ ) levels in biological samples. The specificity of the  $PGE_2$  ELISA was investigated using authentic  $PGE_2$  and fatty acids which, based on their structure, might be anticipated to compete with  $PGE_2$  for binding to antibodies against  $PGE_2$ . The anti-  $PGE_2$  antibody showed virtually no cross-reactivity with other tested eicosanoids (see Table 1).

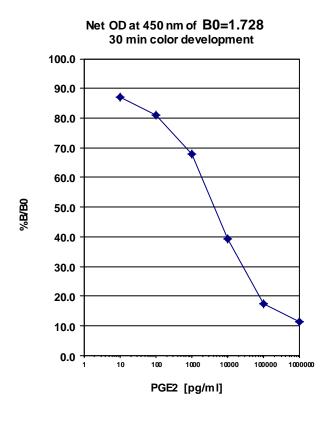

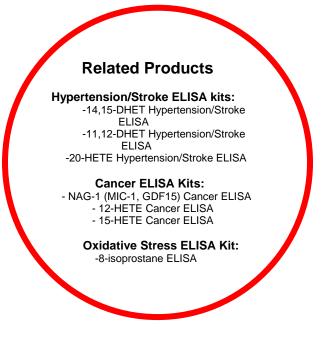


Table1. Specificity of PGE<sub>2</sub> ELISA

Eicosanoids	% Binding of control
PGE <sub>2</sub>	100.00
20-hydroxy PGE <sub>2</sub>	1.16
11-deoxy PGE <sub>2</sub>	<0.01
20-HETÉ	<0.01
12-HETE	<0.01
15-HETE	<0.01
8-isoprostane	<0.01
14,15-DHET	<0.01
11,12-DHET	<0.01